



medical **ECONET**
GERMANY



Only 3,6 kg

Wireless hand switch

Focal spot of 0,8 mm



mex+ AirRay mini

Compact and ultra lightweight portable X-ray generator

- **HYBRID** – Device can be operated by internal battery or by external power supply
- Only 3,6 kg for practical handling in mobile usage
- Equipped with high-performance battery (Lithium-Ion)
- Up to 118 exposures with only one battery charge
- Charging without removing the battery
- Digital Distance Measurement
- 9 programmable pre-set technique slots for kV and mAs
- Built-in ultra capacitor for steady power output
- Tripod mounting thread (1/4" - 20 UNC)

Specifications

Modell	meX+AirRay mini
Output	
Tube Vol. / Current	90kV / 10 mA
Voltage range	50 - 90 kV in 1 kV steps
mAs Range	0,1 mAs -10 mAs 19 steps
Max. output	500 W at 90 kV
X-ray Tube	
Application	OX / 80 CEI
Focal Spot	0,8 mm
Collimator	
Max. Light field: 46 cm x 46 cm @ 100 cm SID	
Lamp	LED
Features	
3,5" LCD Touch display, Digital Distance Measurement	
Exposures full charge	max. 118 at 60 kV / 10 mAs
PROM Memory	9 APR
Weight	3,6 kg
Dimensions	249 x 194 x 194 mm

Standard Accessory

Hand switch, wireless hand switch, 1x rechargeable battery, Carrying Bag



Display Features

Favourites button:
Dose value

Rotate the display



Digital Distance Display


Special feature

Tripod mounting thread (1/4" - 20 UNC)



Enables fast and easy fixation at standard camera tripods.

Dimensions Generator (mm)

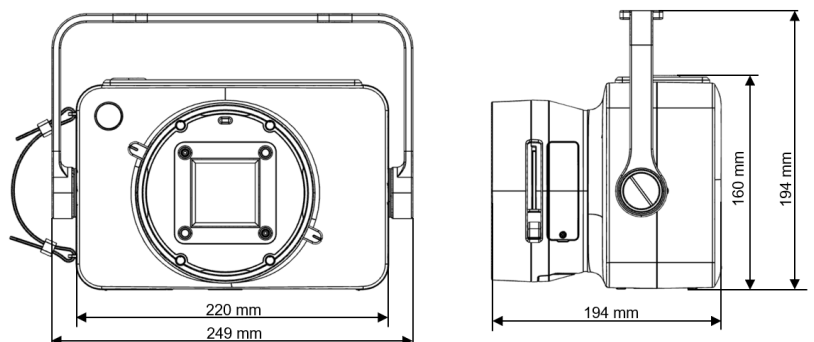

Optional: Hand strap for safe handling.

CE 0123

medical **ECONET**

Im Erlengrund 16, 46149 Oberhausen - Germany

t: +49 (0) 208 377 890-0 info@medical-econet.com
f: +49 (0) 208 377 890-55 www.medical-econet.com
facebook.com/medicaleconetgermany

Distributor

*Subject to change